


Calipers Standard, Extended and Head/Neck/ Mammo Models


Huestis Medical calipers are designed to be the most accurate calipers in radiology. Their detailed construction is precisely controlled so that scales can be calibrated and verified in 1 mm increments to help prevent measurement miscalculations.

Interrelated parts help prevent looseness at the connecting joints to lessen the chance of miscalculation or excess radiation. Precision machined and gold anodized, Huestis Medical calipers offer reliable performance and lasting value.

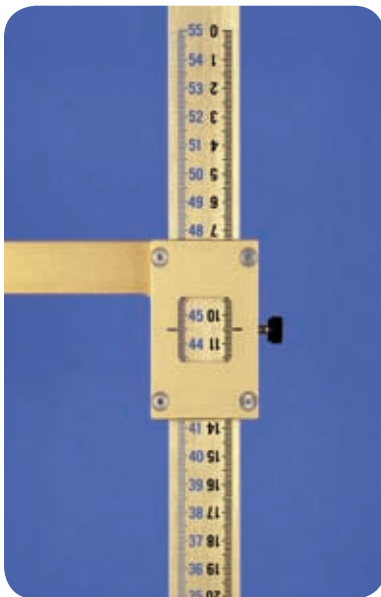

Detailed Construction

Accurately calibrated scales have reversed orientation color-coded numbers for correct reading of various measurements. Recessed indicators protect scales from wear.

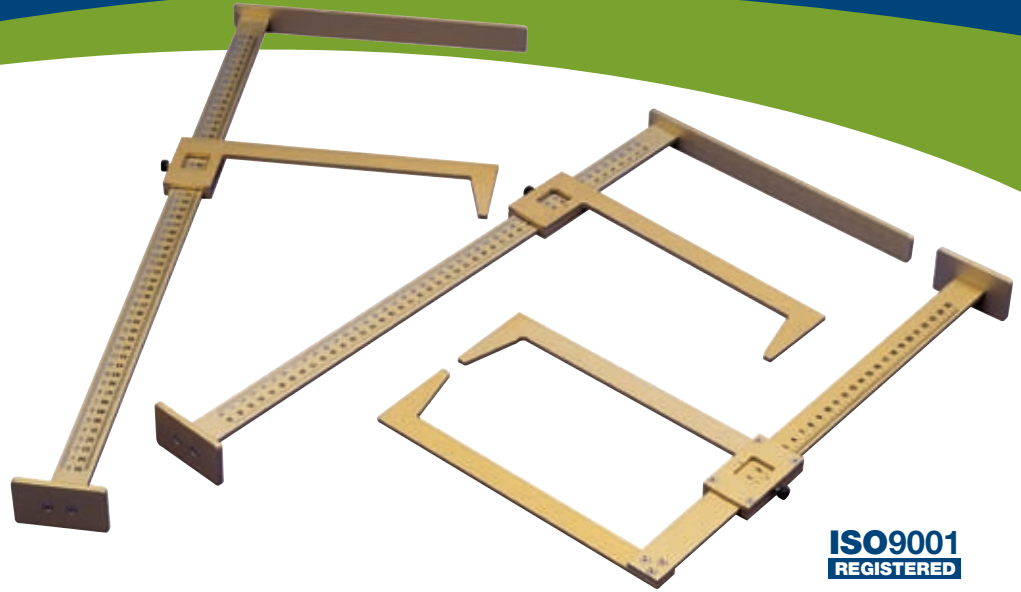

**ISO9001
REGISTERED**

Features:

Accurate Design

- Detailed engineering utilizes minimal fastener points to ensure rigidity.
- Scales are accurately calibrated to 1 mm increments.

Easy to Read

- Reversed orientation scales help prevent reading the wrong scale or transposing numbers on the scale.
- Off-table, point-to-point scales read from the top down.
- Tabletop scales read from the table up, in a different color.
- Recessed indicators protect scales from premature wear.

Versatile Measurements

- Central-axis
- Off-axis
- Laterals
- Obliques
- Off-table

Ease Of Use

- Wide base ensures perpendicularity to table top.
- Specially shaped pointer(s) reach around extended body contours.
- Easily adjusted thumb screw grips firmly at any position on scale.

Precision Machined

- Precision machined, lightweight aluminum ensures accurate tolerances and lasting durability.
- Interrelated parts help prevent looseness at connecting joints.

Hand Assembled

- Hand-tapped, counterbored, and assembled with steel screws, to provide lasting strength.

Gold Anodized

- Protective gold anodizing process prevents deterioration over time.

Standard Body Calipers

- Width: 25.4 cm (10 in)
- Length: 52.8 cm (20.8 in)
- Measures: 0 – 45 cm (0 – 17.7 in)

Extended Body Calipers

- Width: 25.4 cm (10 in)
- Length: 62.3 cm (24.5 in)
- Measures: 0 – 55 cm (0 – 21.7 in)

Head/ Neck/ Mammo

Body Calipers

- Width: 24.7 cm (9.75 in)
- Length: 39.1 cm (15.4 in)
- Measures: 0 – 25 cm (0 – 9.8 in)